



Seahorse Plus Toileting & Shower Chair



The Seahorse Plus enables comfortable and safe toileting and bathing for children aged 5+ to adult

Designed & Developed by NRS Healthcare



The Seahorse Plus inspires confidence and makes toileting and bathing a more manageable experience for both the carer and user.

Suitable for ambulant and non-ambulant users who require postural support, as well as users with mild to severe cognitive impairment or reduced motor function.

Fully adjustable for user security and comfort the Seahorse Plus can quickly and easily be positioned for use over a toilet or bath, as well as in the shower. It can also be used with a suitable hoist. Please contact the Customer Service team for advice on compatible hoists.

Available in three sizes: small, medium and large, with a range of comfortable PU liners to accommodate children as they grow, from 5 years to adult.

A range of optional accessories are available, including; liners, harnesses and anterior supports.

Product features



Tilt-in-Space: The innovative gas-assisted tilting mechanism enables 40° of rotation. To recline the seat simply press the red handle buttons then guide the user from a sitting position to fully reclined for safe, comfortable showering.



100mm seat height adjustment range enables optimal set up to suit most toilets and carer heights.







Comfortable removable armrests that can be substituted for the optional flip-up anterior supports.



Height and depth adjustable headrest with a lateral adjustment feature, to allow for variations in user head position.



Height adjustable footrests with soft heel guard and foot straps that feature a new lockable, flip up design to make entry and exit easy.



Multi adjustable lap-strap with central buckle and tightening feature.



Removable pommel / splash-guard aid user security and positioning.



Removable commode pan

Seahorse Plus Optional Accessories

A range of optional accessories have been designed around the user and carer to offer additional safety, security and maximum comfort.



Seahorse Plus Liners:

Easy to fit, these soft PU foam liners have been designed for comfort and hygiene. The liners can be supplied in different thicknesses, allowing the chair to accommodate the user as they grow.

Small Thick Liner: SH013
Small Thin Liner: SH017
Medium Thick Liner: SH014
Medium Thin Liner: SH018
Large Extra Thick Liner: SH016
Large Thick Liner: SH015
Large Thin Liner: SH019



Chest Harness:

Optional Poziform[™] Chest Harness for users requiring more secure positioning when using the Seahorse Plus Toileting and Shower Chair.

Small: SH020 Medium: SH021 Large: SH022

Extra Large: SH023



Butterfly Harness:

Particularly suitable for women the optional Poziform™ Butterfly Harness offers more secure positioning when using the Seahorse Plus Toileting and Shower Chair.

Small: SH024 Medium: SH025 Large: SH026

Extra Large: SH027







The Anterior Support helps with positioning and aids the user during the toileting process.

Small: SH044


Medium / Large: SH045









Seahorse Plus Dimensions

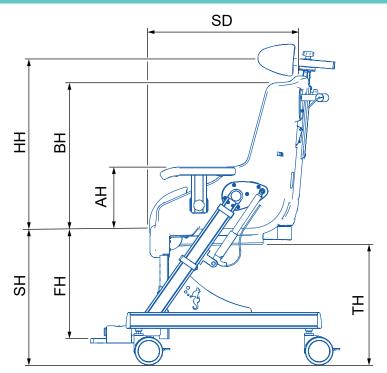

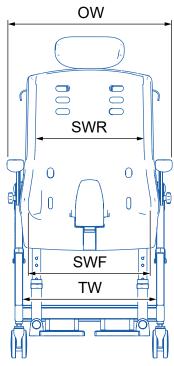

		Small Seahorse Plus SH001						Medium		
Feature	Ref	Thick Liner		Thin Liner		No		Thick Liner		
. 55.41.5		SH013		SH017		Liner		SH014		
Seat Depth	SD	327mm	12¾"	355mm	14"	374mm	14³/₄"	425mm	16³/₄"	
Rear Width of Seat	SWR	258mm	10"	294mm	111/2"	332mm	13"	303mm	12"	
Front Width of Seat	SWF	294mm	11½"	334mm	13¼"	365mm	141/4"	353mm	14"	
Height of Back	ВН	391mm	15½"	391mm	15½"	414mm	16¼"	490mm	191/4"	
Min Headrest height	НН	456mm	18"	456mm	18"	475mm	18½"	562mm	22"	
Max Headrest height	НН	631mm	25"	631mm	25"	650mm	25½"	737mm	29"	
Min height of Seat	SH	540mm	21½"	540mm	21½"	521mm	20½"	521mm	20½"	
Max height of Seat	SH	648mm	25½"	648mm	25½"	629mm	24¾"	625mm	241/2"	
Min height of Footrests	FH	275mm	10¾"	275mm	10¾"	255mm	10"	325mm	12¾"	
Max height of Footrests	FH	375mm	14¾"	375mm	14¾"	355mm	14"	425mm	16¾"	
Total width	OW	532mm	21"	532mm	21"	532mm	21"	586mm	23"	
Total depth (upright)	OD	807mm	31¾"	807mm	31¾"	807mm	31¾"	860mm	33¾"	
Max turning circle (tilted)	TC	1209mm	47½"	1209mm	471/2"	1209mm	471/2"	1413mm	55¾"	
Min height over toilet	TH	461mm	18"	461mm	18"	461mm	18"	442mm	17½"	
Max height over toilet	TH	569mm	22½"	569mm	22½"	569mm	22½"	551mm	21¾"	
Min width over toilet	TW	426mm	16¾"	426mm	16¾"	426mm	16¾"	478mm	18¾"	
Armrest to Seat	АН	179mm	7"	179mm	7"	199mm	8"	214mm	81/2"	
Seat angle	SA	0°-40°		0°-40°		0°-4	0°-40°		0°-40°	
Max User Weight		64kg	10st	64kg	10st	64kg	10st	102kg	16st	
Weight (including liners)		21kg	46.3lbs	20.2kg	44.5lbs	18.7kg	41.2lbs	26.1kg	57.5lbs	



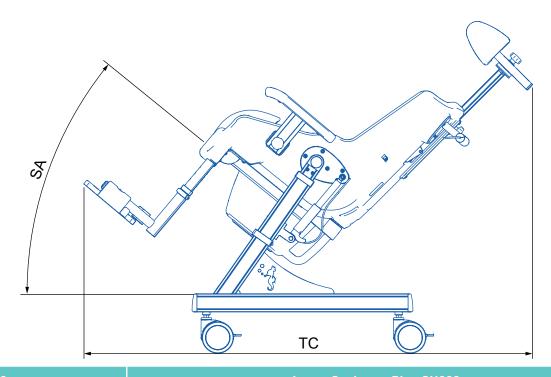

Seahorse Plus SH002			Large Seahorse Plus SH003								
Thin Liner No		Extra Thick Liner		Thick Liner		Thin Liner		No			
SH018 L		Lin	er	SH016		SH015		SH019		Liner	
457mm	18"	470mm	18½"	455mm	18"	455mm	18"	455mm	18"	488mm	191/4"
339mm	13½"	378mm	15"	375mm	14¾"	411mm	16"	447mm	17½"	477mm	19"
393mm	15½"	423mm	16¾"	420mm	16½"	460mm	18"	500mm	19¾"	530mm	20¾"
490mm	191⁄4"	514mm	201/4"	549mm	21½"	549mm	21½"	549mm	21½"	572mm	22½"
562mm	22"	584mm	23"	620mm	24½"	620mm	241/2"	620mm	24½"	641mm	25"
737mm	29"	759mm	30"	795mm	31½"	795mm	31½"	795mm	31½"	816mm	32"
521mm	20½"	501mm	191/2"	521mm	20½"	521mm	20½"	521mm	20½"	501mm	19½"
625mm	241/2"	609mm	24"	629mm	24¾"	629mm	243/4"	629mm	24¾"	609mm	24"
325mm	12¾"	305mm	12"	325mm	12¾"	325mm	12¾"	325mm	12¾"	305mm	12"
425mm	16¾"	405mm	16"	450mm	17¾"	450mm	17¾"	450mm	17¾"	430mm	17"
586mm	23"	586mm	23"	692mm	27"	692mm	27"	692mm	27"	692mm	27"
860mm	33¾"	860mm	33¾"	877mm	34½"	877mm	34½"	877mm	34½"	877mm	34½"
1413mm	55¾"	1413mm	55¾"	1488mm	58½"	1488mm	58½"	1488mm	58½"	1488mm	58½"
442mm	171/2"	442mm	171/2"	442mm	171/2"	442mm	171/2"	442mm	171/2"	442mm	171/2"
551mm	21¾"	551mm	21¾"	550mm	21¾"	550mm	21¾"	550mm	21¾"	550mm	21¾"
478mm	18¾"	478mm	18¾"	592mm	231/4"	592mm	231/4"	592mm	231/4"	592mm	23¼"
214mm	81/2"	234mm	9"	214mm	81/2"	214mm	81/2"	214mm	81/2"	234mm	9"
0°-40° 0°-40°		0°-40°		0°-40°		0°-40°		0°-40°			
102kg	16st	102kg	16st	152kg	24st	152kg	24st	152kg	24st	152kg	24st
24.7kg	54.5lbs	22.3kg	49.2lbs	29.0kg	63.9lbs	28.3kg	62.4lbs	27.7kg	61.1lbs	24.8kg	54.7lbs

Seahorse Plus Details



The Seahorse Plus has been created to support boys and girls as they grow. The chest harness is designed for boys and the butterfly harness for girls, for maximum comfort as they develop.

Seahorse Plus Toileting & Shower Chair

Small (SH001) Medium (SH002) Large (SH004)

Standard Features:

Seahorse Plus Toileting & Shower Chair

Height and depth adjustable headrest

Height adjustable flip-up footrests

Lap-strap

Removable pommel

Removable armrests

Removable commode pan

Accessories								
	Small	Medium	Large	Extra Large				
	Code	Code	Code	Code				
Thick liners	SH013	SH014	SH015	SH016				
Thin liners	SH017	SH018	SH019					
Chest harness	SH020	SH021	SH022	SH023				
Chest harness (no buckle)	SH028	SH029	SH030					
Butterfly harness	SH024	SH025	SH026	SH027				
Anterior supports	SH044	SH045	SH045					

For more information or to place an order, please contact the Smirthwaite Customer Service Team:

Telephone: 01626 835552 Email: info@smirthwaite.co.uk

www.smirthwaite.co.uk

