


# Qai

## paediatric car seat

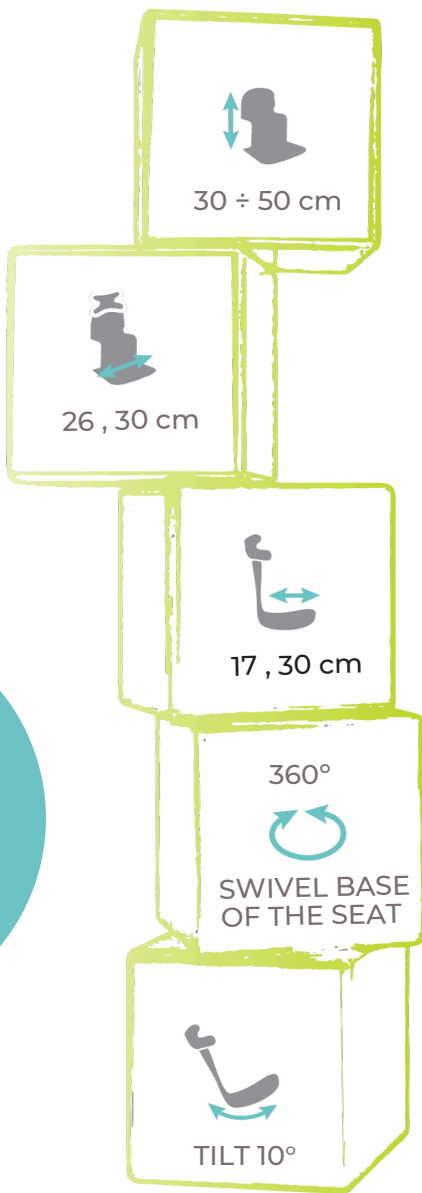

simon

i-Size  
Universal ISOFIX  
61-105cm/≤18kg  
i-Size  
Booster Seat  
100cm-150cm



040365 S  
Regulation No.129/04

**All children have the right to travel safely and comfortably**

Introducing the QAI paediatric car seat; meticulously designed, tested, certified, and approved in compliance with the European standard for child restraint systems used in motor vehicles (ECER129/04 S).

The QAI has been authorized for use as a child restraint system for special needs designated as "S." We have rigorously tested and certified all additional accessories, including lateral supports, hip pads, headrests, and footrests, ensuring compliance with stringent safety regulations. Over 40 crash tests have been conducted to meet the highest safety standards outlined in the ECER129/04 S regulations. We prioritise safety without compromise, as your child deserves the best available solutions.

QAI is suitable for children 61-150 cm of height.



### Important features



The seat can be rotated 360° without any restrictions. The function secures easy access to get into and out of the car.



Rearward facing position for children 61-105cm of height. Forward facing position is allowed from 76cm.



Wide range of positioning and support systems: rigid lateral supports, harness, headrest and abduction wedge.



Reclining positions (tilt – around 10°) and adjustable height of the backrest.



Inner cushions and additional hip pads to reduce the depth and width of the seat. Designed for small and slim kids.



Footrest platform with height and angle adjustment.

### Technical data

Backrest height	30 ÷ 50 cm
Headrest height	52 ÷ 73 cm
Seat depth	30 cm / with cushions 17 cm
Seat depth	30 cm / with cushions 26 cm
Footrest platform width and depth	30 x 19 cm

### Colours

